

initial Zirconia Disk Multilayer Elite

The Top Choice for the Highest Standards



Since 1921
100 years of Quality in Dental



The top choice for the highest standards

Engineered to deliver the perfect balance of high strength and exceptional aesthetics, GC Initial® Zirconia Disk Multilayer Elite is the ideal choice for both anterior and posterior restorations.




SUPERIOR AESTHETICS

True-to-nature shade and translucency gradient for anterior and posterior restorations.



HIGH STRENGTH

Classified as a Class V according to DIN EN ISO 6872 with a mean value of > 1.100 MPa.

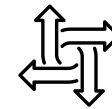

OMNI-FUNCTIONAL

Versatile and suitable for a broad range of clinical indications. Fast sintering or regular sintering.


Natural gradient

Stunningly natural: a seamless gradient from cervical to incisal, in shade as well as translucency gradients, mimicking the natural tooth structure in the best possible way.

Translucency >49% Strength 800 MPa		5Y INCISAL 20%
		4-5Y INTERLAYER 1 15%
		4-5Y INTERLAYER 2 15%
		4-5Y INTERLAYER 3 15%
Translucency >45% Strength 1.150 MPa		4Y CERVICAL 35%

Physical properties

- Zirconium dioxide Type II Class V DIN EN ISO 6871
- Flexural strength: typical mean value ≥ 1.100 MPa
- High translucency gradient from >45% to >49%
- Fast sintering for crowns and bridges up to 3 units

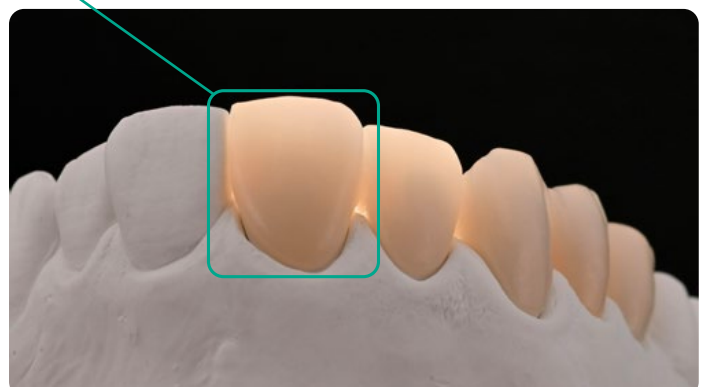


Made for qualitative productivity

Suitable for high productivity without compromising on quality.

SCAN & DESIGN

Gesture-controlled scanner with high precision and smart scan features

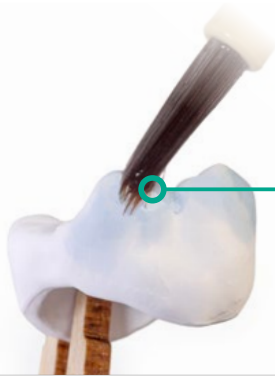

OPTIONAL

INNER STAIN with Coloring Liquid

Refine with effect colors for the ultimate sophistication.

OPTIONAL

MICROLAYERING with GC Initial® IQ ONE SQIN

3D PAINT with GC Initial® IQ Lustre Pastes ONE

Color and fluorescence - either for purely monolithic work or with micro-cutback in combination with GC Initial® IQ ONE SQIN.



FORM & TEXTURE with GC Initial® IQ ONE SQIN

Vitality, natural texture and gloss within a microlayer.



OPTIONAL

CLASSIC VENEERING with GC Initial® Zr-FS


Unsurpassed optical effects and high stability after multiple firings for frameworks with classic cutback.



OR

CEMENT with GC FujiCEM® Evolve

The preferred choice for zirconia, easy to use and ideal for juxtagingival or subgingival margins.

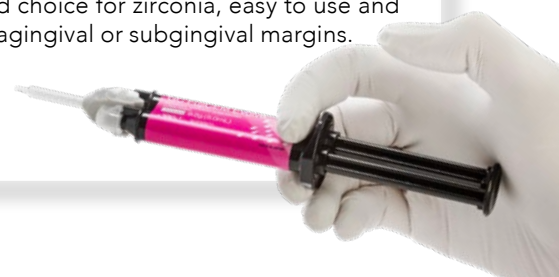

initial Zirconia Disk Multilayer Elite

The Top Choice for the Highest Standards

Benefits at a glance

- Superior aesthetics with a combined shade and translucency gradient.
- High strength throughout all layers, making it suitable for a wide range of indications.
- Limits the required finishing with a framework in the right color tone.
- Optimal synergy with GC Initial® Zr-FS and the GC Initial® IQ ONE SQIN concept for high-end individualization.
- All 16 V-shades and 3 bleach types are available in 14, 16, 20 & 25 mm thick disks with a diameter of 98.5 mm.



Perfect synergies



GC Initial® IQ ONE SQIN



GC Initial® Zr-FS



GC Initial® IQ Lustre Pastes ONE

GC America Inc.

3737 West 127th Street
Alsip, IL 60803
800.323.7063
www.gc.dental/america/
© 2025 GC America Inc.