



Two modes, Two devices
Now in One



Labolight DUO from GC

LED Dual Mode Light Curing Unit
for indirect composite techniques



GC Labolight DUO

The days when you needed two curing devices in your laboratory (one for intermediate and another for final curing) are now a thing of the past. GC is now offering a state of the art multi-functional light curing device for the technician combining 2 curing modes: pre-curing (step mode) and final curing (full mode)

Latest LED technology inside


Equipped with double wave length LED technology the Labolight DUO can be used to cure any of GC's composites in a secure and durable way. It offers a wavelength range of 380nm – 510nm with spectrum ranges peaks of:

- 465 nm-485 nm (12 Blue LEDs)
- 390 nm-400 nm (3 Violet LEDs)

This technology ensures optimal hardening of all light cured dental materials, while the high power outlet reduces light curing cycles.

Automated rotary system

The reflective plate distributes all light effectively and exposes your works from all sides. The curing stand carefully positions the objects during all lightcuring cycles.

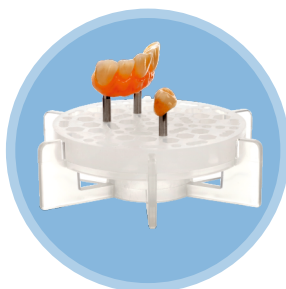

Compact & ergonomic design

The GC Labolight DUO is not only capable of handling most of your work, it looks good as well. The award winning sleek, contemporary styling with smooth surfaces allows for easy cleaning and reduces the space it takes in your lab.


This design serves the usability: the interface is very simple and intuitive and the inner parts are easily accessible thanks to the wide opening.



Step mode



Full mode

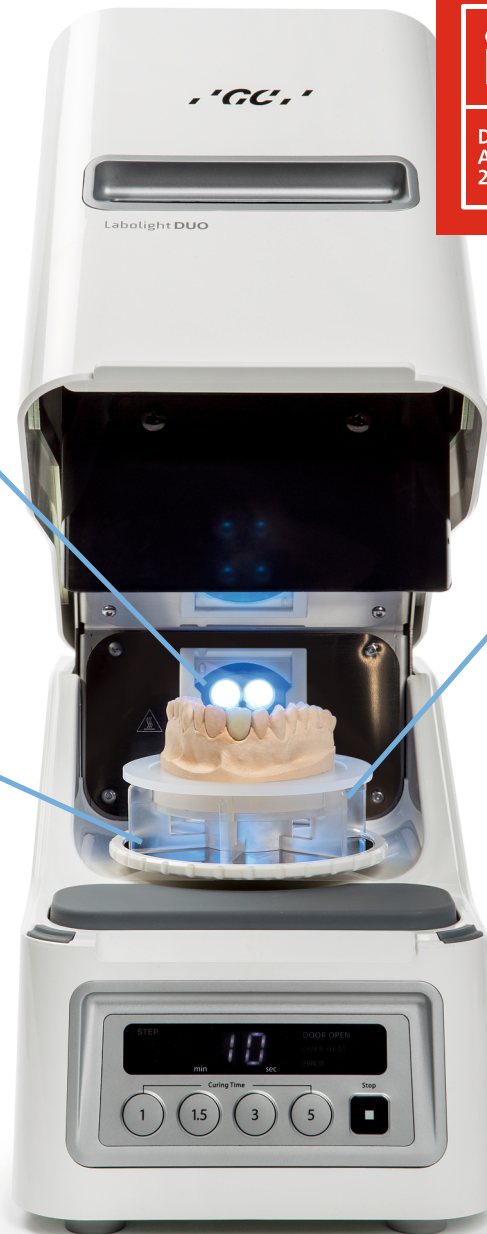

GC Associates:

NORTH INDIA: Amit Sharma - 9818956660, Pramod Sati - 9592514153, **WEST INDIA:** Ganesh Shinde - 9833233172, **EAST INDIA:** Pradiip Sengupta - 9038000265/9830079266, **AP TELANGANA:** Md. Shareef - 9848450321, **KARNATAKA:** Srikanth - 8722272222

NATIONAL SALES MANAGER: Harisankar Radhakrishnan - 9100932652



GC India Dental Pvt. Ltd.

Plot No. 233, Phase III, IDA Pashamylaram,
Pattancheru Mandal, Medak District - TS 502 307
Tel. : +91 8455 224 844 Toll Free : 1800-425-3132

Email : info@gcindiadental.com

Website : www.gcindiadental.com

Connect us on:

- <https://www.facebook.com/GcIndiaDental1>
- https://twitter.com/GCIndia_Dental
- [in.linkedin.com/pub/gc-india-dental/a2/435/63/](https://www.linkedin.com/pub/gc-india-dental/a2/435/63/)
- <http://www.pinterest.com/gcindiadental/>
- <http://plus.google.com/+gcindiadental>

GC INDIA App



for Android



for iPhone