



## Aligner Accessories



Global Experts in Orthodontics

03

**GC Orthodontics**  
GLOBAL EXPERTS IN ORTHODONTICS

## Aligner Accessories

Aligner Pliers	129
G-ænial™ Universal Flo	130
G-ænial™ Universal Injectable	130
G-Premio BOND™	130
MI Paste® ONE Toothpaste	131
MI Paste® & MI Paste Plus®	131
Aligner Chewies	132
Retainer Brite	132
Retainer Cases	132
Bondable Lingual Loops/Surgical Buttons	132
DentaSonic Safe Interproximal Reduction System	132
GC Metal Finishing Strips/Interproximal Strips	132

# Aligner Pliers




**NEW!**

## Hole Punch Pliers

The Hole Punch Plier creates half-moon cutouts to facilitate the use of bondable buttons for elastics with aligners. It eliminates the need to attach buttons to plastic trays. It is also helpful for trimming aligners to reduce impingement and for patient comfort.


		Length	PART NUMBER
		125mm	GCOA-731

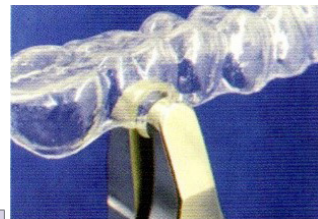

## Horizontal Indent Pliers



The Horizontal Indent Plier is used to create an indentation to accent individual root torque. It can also be used to place an indent into a specific tooth's undercut or below a bonded attachment to enhance retention of clear aligners and retainers.

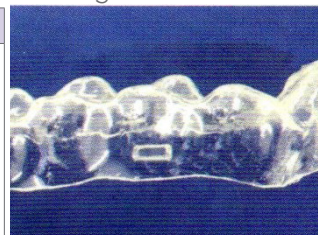
		Length	PART NUMBER
		118mm	GCOA-732



## Vertical Indent Pliers

The Vertical Indent Plier is used to create an indentation at the mesial or distal line angles on the lingual and/or buccal of a specific tooth to emphasize over-correction of rotations. Vertical indents on the lingual and buccal surfaces on opposing sides of a tooth can be created. In/out discrepancies can also be addressed with this plier.

		Length	PART NUMBER
		118mm	GCOA-733

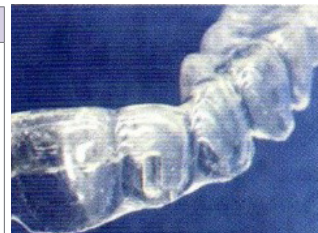

## Teardrop Hole Pliers


The Teardrop Hole Plier is used to create a notch at the gingival margin of clear aligners where elastic hooks are needed. It provides a reservoir to hold elastics in place as the aligner is seated, which helps the patient with placement of elastics.

		Length	PART NUMBER
		118mm	GCOA-734




## G-aenial™ Universal Flo

G-aenial Universal Flo was released in 2011 as the first universal application, high strength flowable composite resin. It has been ideal for creating aligner attachments with excellent clinical outcomes (exceptional aesthetics, no flash, ideal flow without being runny, high strength and no porosities) and most importantly a significant reduction in application time helping further improve the cost effectiveness of the aligner treatment procedure.

	Description	Quantity	PART NUMBER
	G-aenial Universal Flo A1	1 Syringe (2 mL/3.4 g per syringe)	004202
	G-aenial Universal Flo A2	1 Syringe (2 mL/3.4 g per syringe)	004203
	Dispensing tips III	30 pcs. (plastic)	004634
	Dispensing tips III, Needle Tip	30 pcs. (metal)	004635



1. Invisalign attachment template  
 2. Teeth are cleaned with pumice, then etched with phosphoric acid for 30 secs, wash and dry.  
 3. Apply G-BOND, leave 5-10 secs, dry with MAXIMUM air for 5 sec's. Light cure for 10 secs.  
 4. G-aenial Universal Flo  
 5. G-aenial Universal Flo, shade AE is used to fill the attachment space in the template. Application is controlled to avoid overfilling.

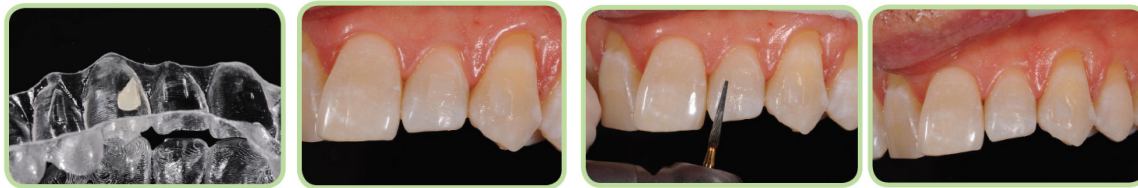


6. Once the template is placed in the mouth, the composite is light cured for 20 seconds, the template is then removed.  
 7. Attachments in position, following template removal.  
 8. Excess flash is removed with a tungsten carbide bur.  
 9. Completed Invisalign attachments after final polishing.

Dr Anthony Mak, Sydney NSW

## G-aenial™ Universal Injectable




G-aenial™ Universal Injectable is a radiopaque, injectable high-strength ultra-fine particle composite. Ideal for aligner attachments, lingual retainers and indirect bonding, it features total control with no slumping, stringing or flash. Resistant to discoloration, it comes in a 2-syringe kit with 10 long dispensing tips and 10 short dispensing tips.

	Description	Quantity	PART NUMBER
	G-aenial Universal Injectable A1	Two 1.7g (1.0mL) syringes, 10 long & short tips	012364
	G-aenial Universal Injectable A2	Two 1.7g (1.0mL) syringes, 10 long & short tips	012365
	G-aenial Universal Injectable JE	Two 1.7g (1.0mL) syringes, 10 long & short tips	012376
	G-aenial Universal Injectable AE	Two 1.7g (1.0mL) syringes, 10 long & short tips	012377
Dispensing tips (Long Needle)	20 pcs.	012281	

## G-Premio Bond™

G-Premio BOND is a universal, 8th generation bonding agent that is compatible with total-etch, self-etch and selective etch techniques providing excellent versatility. Fast application (10 seconds). Stability with extended working times (up to 5 minutes).

	Description	Quantity	PART NUMBER
	G-Premio Bond Kit	1 bottle (5mL), 50 disposable applicators, 20 disposable dispensing dishes	009275
	G-Premio Bond Refill	1 x 5 mL bottle	009277
	G-Premio Bond Unit Dose Refill	50 unit doses (0.1 mL) 50 disposable applicators (fine)	009276

# MI Paste




## MI Paste® ONE Toothpaste

MI Paste® ONE is the newest member of the MI Paste® family as a 2-in1 application (cleaning and treatment) containing RECALDENT™ (CPP-ACP) and fluoride, in a cool mint flavor!

MI Paste® ONE is used for your patients to safely:

1. Remineralize
2. Rebuild Enamel
3. Relieve Sensitivity.



	Description	Quantity	PART NUMBER
	MI Paste One Toothpaste (Mint)	One case (ten tubes per case).	437000

## MI Paste® & MI Paste Plus®

Topical Tooth Crème Containing Calcium, Phosphate & Fluoride




MI Paste and MI Paste Plus are the only products for professional use containing the active ingredient RECALDENT™ (CPP-ACP), a special milk-derived protein that has a unique ability to release bio-available calcium and phosphate (and fluoride in MI Paste Plus) to tooth surfaces. When it comes to maintaining a healthy smile, these key ingredients can make a big difference in overall health.

MI Paste Plus offers the same benefits of MI Paste, but is enhanced with a patented form of fluoride (900 ppm). MI Paste and MI Paste Plus are safe and easy-to-use both in office and at-home. They are both available in five delicious flavors: melon, mint, strawberry, tutti-frutti and vanilla.


### What is the difference between MI Paste & MI Paste Plus?

MI Paste does not contain fluoride and is recommended for children under six years of age and pregnant women. MI Paste Plus offers the same great benefits of regular MI Paste; it is enhanced with a patented form of fluoride.

### MI Paste

	Description	Quantity	PART NUMBER
	MI Paste (Assorted Package)	Contains 10 tubes, 2 of each flavor: melon, mint, strawberry, tutti frutti and vanilla (40 g; 35 mL each)	422265
	MI Paste (Mint)	Contains: 10 tubes (40 g; 35 mL each)	423679
	MI Paste (Strawberry)	Contains: 10 tubes (40 g; 35 mL each)	424505


### MI Paste Plus

	Description	Quantity	PART NUMBER
	MI Paste Plus (Assorted Pkg)	Contains 10 tubes, 2 of each flavor: melon, mint, strawberry, tutti frutti and vanilla (40 g; 35 mL each)	422614
	MI Paste Plus (Mint)	Contains: 10 tubes (40 g; 35 mL each)	422621
	MI Paste Plus (Strawberry)	Contains: 10 tubes (40 g; 35 mL each)	422886
	MI Paste Plus (Vanilla)	Contains: 10 tubes (40 g; 35 mL each)	422888


CPP-ACP was developed at the School of Dental Science at the University of Melbourne Victoria/Australia. RECALDENT™ is a trade mark used under license. RECALDENT™ (CPP-ACP) is derived from milk casein. Do not use on patients with a milk protein and/or hydroxybenzoates allergy.

# Aligner Accessories


## Ligature and Band Scissor (Ideal for trimming invisible retainer material, cutting bands, straight band material and fine wires)

	Description	Quantity	PART NUMBER
	Ligature and Band Scissor, Light - Straight	each	BECK500-203
	Ligature and Band Scissor, Light - Curved	each	BECK500-207


## Aligner Chewies, Retainer Brite, and Retainer Cases

	Description	Quantity	PART NUMBER
	Aligner Chewies (Unscented)	20/pkg	55-0160-0000
	Retainer Brite (3 month supply)	96/pkg	55-0212-0000
	Retainer Cases (15 Assorted Colors) Dimensions: 1" deep, 2-1/2" wide, 3" front to back	50/pkg	55-0210-0035


## Seagull Buttons

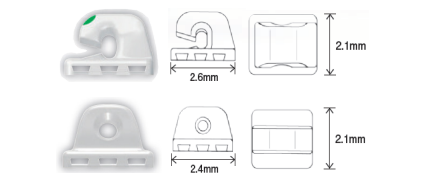	Description	Quantity	PART NUMBER
	Bondable Seagull Button, Open	20/pkg	44-2130-0000
	Bondable Seagull Button, Closed	20/pkg	44-2140-0000


## Bondable Lingual Loops/Surgical Buttons

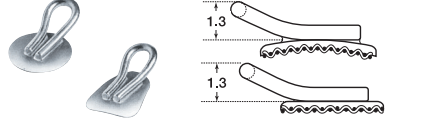	Description	Quantity	PART NUMBER
	Lingual Loop For Curved Teeth	10/pkg	42-0030-0000
	Lingual Loop For Flat Teeth	10/pkg	42-0040-0000

## DentaSonic Safe Interproximal Reduction System (see pages 165-166 for full line and options)

	Description	Quantity	PART NUMBER
	<b>DentaSonic Profi Kit 2 - 4:1 Water Cooled Hand Piece</b> Contains water cooled hand piece, additional head, 4 single and 4 double sided strips, polishing strips, debond polishers, measuring gauge, extraction tip and spray nozzle	Kit	55-BD02-0000
	<b>DentaSonic Basic Set - 4:1 Non Water Cooled HP</b> Contains non-water cooled hand piece, measuring gauge, 4 assorted strips, extraction tip and spray nozzle	Kit	55-BD04-0000

## GC Metal Finishing Strips/Interproximal Strips (Contains 12 strips with the same Mesh and Width)

	Description	Quantity	PART NUMBER	Description	Quantity	PART NUMBER
	Extra Coarse; Red; 2.6 mm	12/pkg	000252	Medium; Green; 2.6 mm	12/pkg	000260
	Extra Coarse; Red; 3.3 mm	12/pkg	000253	Medium; Green; 3.3 mm	12/pkg	000261
	Extra Coarse; Red; 4.0 mm	12/pkg	000254	Medium; Green; 4.0 mm	12/pkg	000262
	Coarse; Blue; 2.6 mm	12/pkg	000256	Fine; Yellow; 2.6 mm	12/pkg	000264
	Coarse; Blue; 3.3 mm	12/pkg	000257	Fine; Yellow; 3.3 mm	12/pkg	000265
	Coarse; Blue; 4.0 mm	12/pkg	000258	Fine; Yellow; 4.0 mm	12/pkg	000266