

Miscellaneous



Global Experts in Orthodontics

03

GC Orthodontics
GLOBAL EXPERTS IN ORTHODONTICS

Miscellaneous

DentaSonic Safe Interproximal Reduction System	165
GC Metal Finishing Strips/Interproximal Strips	167
Cheek Retractors	167
Photographic Mirrors	167
Microbrushes	167
Retainer Cases	168
Swivel Safety Keys (RPE)	168
Patient Wax and Silicone	168
Elastibites	168
Aligner Chewies	168
Retainer Brite	168
Bracket Cards	169
Elastics Placers	169
Aluminum Oxide	169
Wire Storage Rack	169
Archwire Markers	169
Silver Solder	169

DentaSonic

Swiss Proximal Restorative Clean-up Tools





The innovative DentaSonic IPR system provides a safe and accurate, minimally invasive enamel reduction system. Eliminating concerns over damaging soft tissue while opening the contact, you now have a system that is comfortable for the patient, and worry-free for orthodontists. The components of the system are autoclaveable for multiple procedures.

The DentaSonic system has several advantages:




1. Predefined thicknesses of the IPR attachments allow you to simply and safely remove exactly the amount of enamel necessary
2. The system is gentle, working with the slow speed hand piece, the oscillating motion has little to no vibration
3. Specially designed contact opening attachments open the contact point quickly, but will not damage a patient's gums
4. A perfect final polish is achieved, with a flexible polishing attachment made of polyester and aluminum oxide for perfect polishing
5. Smooth and quiet contra-angle hand piece with water cooling provides a gentle, oscillating movement at a stroke of 1.2mm

Starter Kits

	Description	Quantity	PART NUMBER
	<p>Starter Set (Profi Kit 2) 4:1 Water Cooling Includes: Water cooling hand piece with 4:1 reduction, 4 DiaStrips - 2 sided, 4 DiaStrips - 1 sided, 2 Polishing Polistrrips - 2 sided, 1 Measuring gauge with 6 thicknesses: 0.1, 0.15, 0.2, 0.3, 0.4, 0.5mm, 1 DentaSonic Twist water cooling hand piece with 4:1 reduction, 1 additional oscillating head, 1 extraction tip, 1 spray nozzle</p>	Starter Set	55-BD02-0000
	<p>Basic Set 4:1 Non Water Cooling Includes: Non water cooling hand piece with 4:1 reduction, 4 assorted Strips, 1 Measuring gauge with 6 thicknesses: 0.1, 0.15, 0.2, 0.3, 0.4, 0.5mm, 1 extraction tip, 1 spray nozzle</p>	Basic Set	55-BD04-0000


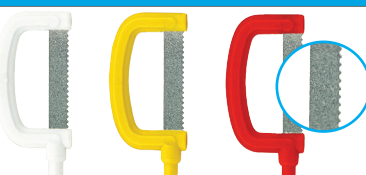
DiaStrip

Oscillating DiaStrips are available 1-sided or 2-sided. Precisely targeted-thickness strips, for controlled IPR from 0.1 to 0.3mm

	Description	Color	Thickness	PART NUMBER	
				5/pkg	12/pkg
	DiaStrip Refill - 1 Side	<input type="checkbox"/> White	0.1 mm, 15 µm	55-FDM0-5115	---
	DiaStrip Refill - 1 Side	<input type="checkbox"/> Yellow	0.15 mm, 25 µm	55-FDM0-5125	---
	DiaStrip Refill - 1 Side	<input type="checkbox"/> Red	0.2 mm, 40 µm	55-FDM0-5140	---
	DiaStrip Refill - 1 Side	<input type="checkbox"/> Blue	0.3 mm, 60 µm	55-FDM0-5160	---
	DiaStrip Refill - 2 Sided	<input type="checkbox"/> White	0.1 mm, 15 µm	55-FDM0-5215	55-FDM0-1215
	DiaStrip Refill - 2 Sided	<input type="checkbox"/> Yellow	0.15 mm, 25 µm	55-FDM0-5225	55-FDM0-1225
	DiaStrip Refill - 2 Sided	<input type="checkbox"/> Red	0.2 mm, 40 µm	55-FDM0-5240	55-FDM0-1240
	DiaStrip Refill - 2 Sided	<input type="checkbox"/> Blue	0.3 mm, 60 µm	55-FDM0-5260	55-FDM0-1260
	DS Assortment - 2 Sided		Assorted	55-FDM0-5200	55-FDM0-1200

DiaSaw

Specially designed for precise contact opening without risk of soft tissue damage



	Description	Color	Thickness	PART NUMBER	
				5/pkg	
	DiaSaw Refill - 2 Sided	<input type="checkbox"/> White	0.1 mm, 15 µm	55-SDM0-5315	
	DiaSaw Refill - 2 Sided	<input type="checkbox"/> Yellow	0.15 mm, 25 µm	55-SDM0-5325	
	DiaSaw Refill - 2 Sided	<input type="checkbox"/> Red	0.2 mm, 40 µm	55-SDM0-5340	
	Assortment - 2 Sided		Assorted	55-SDM0-5300	

DentaSonic

Swiss Proximal Restorative Clean-up Tools



DiaStrip Aligner Set

Comprehensive set of strips most commonly needed in aligner treatment (assortment of 5)

 DiaStrip	Description	Quantity	PART NUMBER
	Includes: 1 .15 DiaSaw, 0.25 and 0.4 DiaStrips, 2 PoliStrips	5/pkg assortment	55-FD0A-5300




PoliStrip

Restorative Polishing Strips that are flexible, made of polyester and aluminum oxide for perfect polishing (package of 12)

 PoliStrip	Description	Quantity	PART NUMBER
	Resin Polistrip for 2-sided polishing	12/pkg	55-FDRA-1215



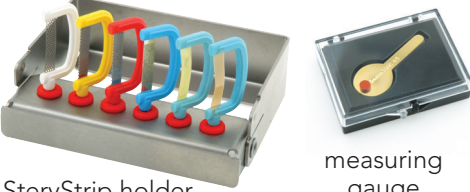

Accessories - SteryStrip Holder and Measuring Gauge

Autoclavable stainless steel tray for cleaning of all strips and attachments (holds 6 strips)/Measuring Gauge for precision IPR

	Description	Quantity	PART NUMBER
	SteryStrip Holder - Holds 6 strips	each	55-TDM0-6000
	Measuring Gauge - 6 thicknesses	each	55-G000-1000







Debonding Polisher

For safe removal of brackets/aligner attachments/adhesives without damage to the tooth surface

	Description	Quantity	PART NUMBER
	Debonding Set (3 point, 3 cup)	6/pkg	55-D000-6B00
	Debonding Set (6 point)	6/pkg	55-D000-6200
	Debonding Set (6 cup)	6/pkg	55-D000-6100

Slow Speed Hand Piece and Additional Parts

4:1 Hand Piece (5,000-7,000rpm), 1.2mm oscillating stroke, autoclavable, to be used with all DentaSonic oscillating attachments


	Description	Quantity	PART NUMBER
	Water Cooling Twist Hand Piece (Fits all air and electric micromotors)	each	55-WD01-0000
	Non-Water Cooling Hand Piece (Fits all air and electric micromotors)	each	55-WD03-0000
	Water Cooling Replacement Head	each	55-KD01-0000
	Non-Water Cooling Replacement Head	each	55-KDR3-0000
	Micro Motor 25'000 rpm	each	55-MD0W-0000







Miscellaneous

Miscellaneous



GC Metal Finishing Strips/Interproximal Strips (Contains 12 single-sided strips with the same Mesh and Width)

	Description	Quantity	PART NUMBER	Description	Quantity	PART NUMBER
	Extra Coarse; Red; 2.6 mm	12/pkg	000252	Medium; Green; 2.6 mm	12/pkg	000260
	Extra Coarse; Red; 3.3 mm	12/pkg	000253	Medium; Green; 3.3 mm	12/pkg	000261
	Extra Coarse; Red; 4.0 mm	12/pkg	000254	Medium; Green; 4.0 mm	12/pkg	000262
	Coarse; Blue; 2.6 mm	12/pkg	000256	Fine; Yellow; 2.6 mm	12/pkg	000264
	Coarse; Blue; 3.3 mm	12/pkg	000257	Fine; Yellow; 3.3 mm	12/pkg	000265
	Coarse; Blue; 4.0 mm	12/pkg	000258	Fine; Yellow; 4.0 mm	12/pkg	000266


Cheek Retractors

	Description	Quantity	PART NUMBER
	Autoclavable Cheek Retractor, Full View, Clear	each	55-0240-0003
	Disposable Cheek Retractor, Medium (Blue)	10/pkg	55-0243-0002
	Disposable Cheek Retractor, Large (Red)	10/pkg	55-0243-0003
	Double-End Large Autoclavable Cheek Retractor, Clear	2/pkg	55-0241-0003
	Photographic Cheek Retractor - Small/Pedo	each	BECK500-322
	Photographic Cheek Retractor - Medium	each	BECK500-323
	Photographic Cheek Retractor - Large	each	BECK500-324

Photographic Mirrors


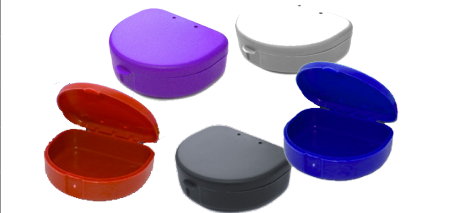
	Description	Quantity	PART NUMBER
	Premium Mirrors (Glass, Set of 4)	4/pkg	55-0170-0000
	Premium Occlusal Mirror (Glass)	each	55-0171-0000

Microbrushes


	Description	Quantity	PART NUMBER
	Microbrushes, Purple	100/pkg	55-0250-0016
	Microbrushes, Green	100/pkg	55-0250-0026
	Microbrushes, Blue	100/pkg	55-0250-0039
	Microbrushes, Peach	100/pkg	55-0250-0046

Miscellaneous



Retainer Cases

	Description	Quantity	PART NUMBER
	Retainer Cases (15 Assorted Colors) Dimensions: 1" deep, 2-1/2" wide, 3" front to back	50/pkg	55-0210-0035
	Retainer Cases - Clear	50/pkg	55-0213-0010
	Retainer Cases - Purple	50/pkg	55-0213-0016
	Retainer Cases - Blue	50/pkg	55-0213-0022
	Retainer Cases - Red	50/pkg	55-0213-0025
	Retainer Cases - Black	50/pkg	55-0213-0027


Swivel Safety Keys

	Description	Quantity	PART NUMBER
	Shilladay Key	10/pkg	55-0145-0000


Patient Wax and Silicone

	Description	Quantity	PART NUMBER
	Patient Relief Wax, Unscented (Red)	50/pkg	55-0120-0000
	Protective Relief Silicone (Orthosil)	50/pkg	55-0121-0000


Elasti-Bites

	Description	Quantity	PART NUMBER
	Elasti-bites (Cherry Flavored)	25/pkg	55-0150-0013

Aligner Chewies




	Description	Quantity	PART NUMBER
	Aligner Chewies (Unscented)	20/pkg	55-0160-0000

Retainer Brite


	Description	Quantity	PART NUMBER
	Retainer Brite (3 month supply)	96/pkg.	55-0212-0000

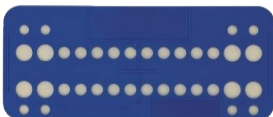


Bracket Cards

	Description	Quantity	PART NUMBER
	Bracket Cards (Sticky)	25/pkg	55-0160-0021
	Single Kit Bracket Cards with Lid	10/pkg	55-7000-0000
	"A-Box" Bracket Setup Cards	10/pkg	55-0300-0000


Elastics Placers

	Description	Quantity	PART NUMBER
	Elastics Placers	100/pkg	55-0140-0021


Aluminum Oxide

	Description	Quantity	PART NUMBER
	Aluminum Oxide, 50 Micron (White)	12 oz bottle	55-0260-0050
	Aluminum Oxide, 90 Micron (Tan)	12 oz bottle	55-0260-0090


Wire Storage Rack

	Description	Quantity	PART NUMBER
	Wire Storage Rack	each	55-7010-0000

Archwire Markers

	Description	Quantity	PART NUMBER
	Archwire Marking Pencils (Red tip)	100/pkg.	55-0125-0025

Silver Solder

	Description	Quantity	PART NUMBER
	Silver Solder 031mm	6 ft.	92-0010-0031